

# MagniLink S Premium 2





# Small, mobile and easy to use

MagniLink S Premium 2 is the ultimate video magnifier for you with high demands on both quality and mobility. It offers superior image quality, modern technology and really smart features. It's also easy to carry, for example, between home and school or a workplace – and can be connected to a PC, Mac, Chromebook, iPad or a monitor. Thus, you can choose the features and options that fit all your needs.





### **Both for distance and reading**

You can easily switch between distance and reading mode. Just turn the camera head and flip up the camera lens and your MagniLink S Premium 2 becomes a distance camera in a second.

### **Easy to handle**

The unit is easy to set up and is ready to use in seconds. Designed with tactile recognition of all main functions, and all parts are made with durability, recyclability and a user friendly feeling in mind.

### **Minimized footprint**

The foldable support legs minimize the footprint, and makes it possible to place the unit tightly alongside a laptop, either on the left or the right side.

### **Superior working height**

The 31 cm working height makes it easy to turn the page on large books and bulky binders.

### **Carefree power supply**

You will never worry about battery strength or misplaced battery chargers as the unit is powered by the computer.



# Make it yours!

All settings are easily controlled by the included software, such as magnification, artificial colors and brightness. This makes it easy to adjust the device according to your needs and preferences.

## **Listen to your texts**

The software included is used to display the image from your MagniLink S on a PC, Mac or Chromebook. To the software you can add TTS (Text To Speech) as an option. With this feature, it is possible to OCR process a text and have it read aloud. The OCR processed text can be displayed in several different display modes, for example line by line, word by word and character by character, while the text is synchronized with speech.

You can OCR process text from an image that you have taken with your MagniLink S, but you can also drag other documents into the program to have them OCR processed. It is also possible to paste text straight into the program to get it read. With text to speech you can increase your reading speed as well as reading endurance.

## **Compatibility**


The software is compatible with other assistive software such as ZoomText and SuperNova.

## **Use with any computer or iPad**

You can use your MagniLink S Premium 2 together Windows, Mac OS, Chrome OS and IOS. If you don't want to connect to a computer you can connect directly to a monitor instead.

## **Split screen functionality**

The software allows you to split the screen image between your video magnifier and your computer – perfect in school- and work environment. This is, for example, perfect in the classroom when the student needs to see what the teacher is doing at the board and at the same time look at the computer.

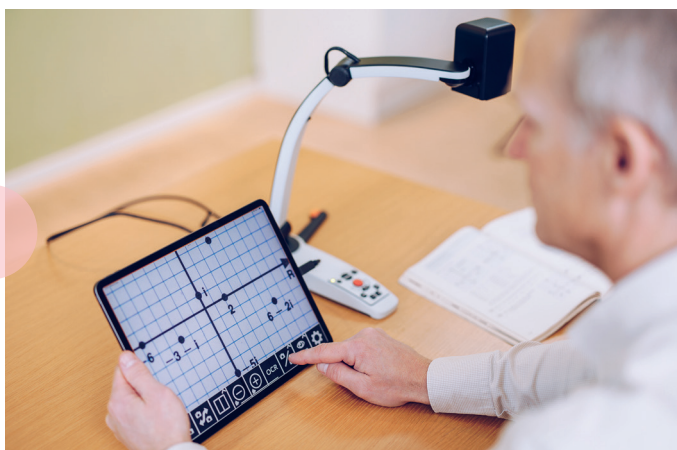





## How to connect to iPad\*



1. Download MagniLink Viewer from App Store.
2. Connect directly to an iPad via USB C.

Recommended to use USB-HUB with minimum 2 ports.

(1 x USB-C PD max. 100W (no data), 1-3 x USB-C 3.2 GEN 2 – max.data speed 10 Gbps / max. 7.5W charging per port)

\* iPad Pro 11" (3rd generation, 2021), iPad Pro 12.9" (5th generation, 2021) and iPad Air (5th generation, 2022) with USB-C-port and IOS 17.0. For questions regarding choice of USB-HUB/AC-adaptor, please contact us (see contact information at last page).



### Smart docking stations

With a docking station you convert your MagniLink S Premium 2 into a stationary system when needed. Dock and undock in seconds to switch between a stationary and a portable system. All docking stations are supplied with an external control panel. The following different docking stations are available:

- **Docking station with X/Y table**

A table in A3 size with friction brakes in both directions.

- **Docking station with table clamp**

Ideal when you have a small work place and want to have the unit at the same spot at the desk.

- **Docking station with plate**

Ideal when you have a small work space but want to be flexible with placement.







Contact us  
and we will answer  
your questions and  
assist with your  
order!

## Technical specifications

### Model

MagniLink S Premium 2

### Connections

USB 3.0 (type C), HDMI

### Camera resolution

HD: 1280 x 720, FHD: 1920 x 1080

### Magnification

FHD: 1.1x-70x\* (23" Monitor)

### Working height

Up to 31 cm (adjustable by the user)

### Colours

Natural colours, 6 artificial colours  
(positive/negative)

### Viewing modes

Read (left/right), distance, mirror

### Compatibility

PC: ZoomText 10.X or later

Dolphin SuperNova 12 or later

### System requirements

**PC:** Windows 7 or later

**Mac:** Mac OS X 10.14 or later

**Chromebook:** Chrome OS 59 or later

**iPad:** iPad Pro 11" (3rd generation, 2021),  
iPad Pro 12.9" (5th generation, 2021) and  
iPad Air (5th generation, 2022) with USB-C-  
port and IOS 17.0.

### Weight

1,9 kg (4.2 lbs)

\* Min. and max. magnification is configurable in the EUC software



# MagniLink S Premium 2

LVI Low Vision International is recognized as one of the world's leading manufacturers of low vision aids with the mission to make daily life easier for people with visual impairments.

We are certified to ISO 14001 & 9001



LVI Low Vision International  
Verkstadsgratan 5  
352 46 Vaxjö, SWEDEN  
Tel: +46 (0)470 72 77 00  
Fax: +46 (0)470 72 77 25  
E-mail: [info@lvi.se](mailto:info@lvi.se)  
[WWW.LVI.SE](http://WWW.LVI.SE)