


Optiva IV

High risk prevention & intensive treatment of wounds up to grade 4.

When bedsores require intensive treatment, only the best will do. The Optiva system offers effective solutions for a wide range of needs in resident/patient management. The mattress allows multiple configurations to offer the best possible pressure redistribution. The system is equipped with an automatic pump, featuring alternating mode, continuous low pressure, pulsation and night mode. Heel cells can be disconnected for 0-pressure in this sensitive area.

Activating the Autofirm Mode provides a stable support surface for nursing procedures or to facilitate transfers in or out of bed. The Lock function improves safety while the audio-visual Alarm guards pressure levels of the system and power supply.

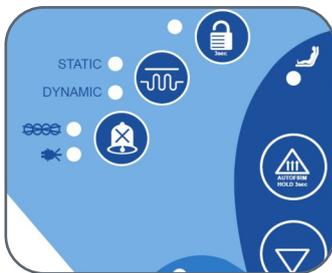

Lock function, Low pressure & Power failure alarm, Autofirm



PU cover with white bottom for inspection of waterproofness



Adjustable hooks



Quick deflation - CPR

<18 sec.*

* Test conditions alternating mode (70-80kg/160-180cm)

Safety Standards

Complies with:
Regulation (EU) 2017 / 745
for medical devices

Protection against Electric shock:
Class II , type BF

Protection of Ingress of Liquids:
IP21

Mode of Operation:
Continuous

Product Specifications

Mattress: Replacement
Width: 85, 90, 100, 120cm
Length: 200, 220cm
Height: 20cm
Patient weight: 40 - 200kg

Pump: 230V / 50Hz
Noise: 0 - 26dB
Cycle: 10 / 20min.

Design, Quality & Service powered

by: **PRIME**
Medical Solutions