

SenSiva


High risk prevention & intensive treatment of wounds up to grade 4.

The ultimate system for pressure redistribution and efficient treatment of bedsores in all care environments. The SenSiva pump detects when a patient is on the mattress and initiates the automatic pressure setting within a few minutes. Once completed the LCD-display will inform the caregiver on patient weight, pressure setting and realtime pressure monitoring of the air cells.

The SenSiva system is both smart and powerful, with adjustable operation modes as Alternating, Continuous Low pressure, Pulsation Therapy and Automatic Fowler pressure boost. Pressure differentiation between cells can be matched with treatment preferences and the system will record usage time for service.




Self sealing connector
for transport

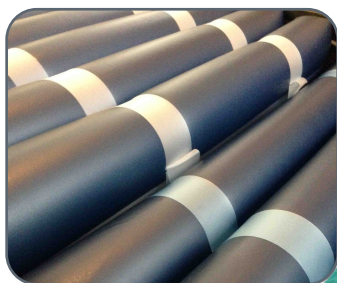

Lower heel cells



Featuring auto weight detection
and pressure differentiation
between cells

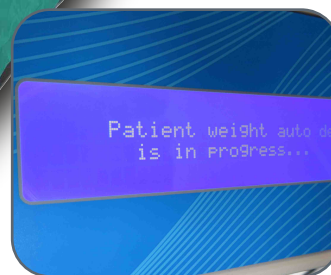

LCD screen and realtime
pressure monitoring

Safety Standards

Complies with:
EC Directive 93/42 Medical Devices
IEC60601-1: 2005 / EN60601-1:2006

Protection against Electric shock:
Class I , type BF

Protection of Ingress of Liquids:
IP21

Mode of Operation:
Continuous

Product Specifications

Mattress: Replacement
Width: 85 or 90cm
Length: 200cm
Height: 20cm
Patient weight: 250kg max.
Weight system: 20kg

Pump: 230V / 50Hz
Noise: 0 - 26dB
Cycle: 10 / 15 / 20min.
Design, Quality & Service powered
by:

PRIME
Medical Solutions