

P4-Adventure P4-Crossover



P4-Adventure
P4-Crossover



4POWER4
Always on track




No path scares it !



Mixed indoor/outdoor use.

Our **P4-Crossover** and **P4-Adventure** wheelchairs share 4 motors and 4-wheel drive (one motor per wheel).

- The **P4-Crossover**, with its fixed rear wheels and power-assisted steering front wheels, offers a turning radius of 1.21m.

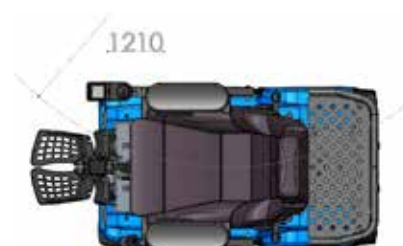

- The **P4-Adventure** also features 4-wheel drive with all-wheel steering controlled by a servo motor, providing one of the best possible turning radius of 0.75m.



Few options

With the **P4**, there are no hidden costs (the basics are included).

They can be customized according to your needs.

Top Case (*)

Practical and clever:

A handle brings it close to you.

It folds away effortlessly with one hand.

Capacity: 39 liters, 3kg load.

Useful in all circumstances

Designed to go through doors, they allow you to move indoors, but more importantly, outdoors in your garden, in the countryside, mountains, seaside, and city.

They will take you easily where other wheelchairs can struggle:

gravel, wet sand (not dry sand), mud, snow, stones, grass, cobblestones... are their playground.



In the countryside



Originally designed for rural environments, with a rust-protected chassis and screws, the **P4** models adapt to many situations such as mud, snow, water (up to 10cm without submerging the motors), stones, etc.

Powerful motorization

4 motors of 350 watts each combined with low-pressure tires provide optimal grip.

For easier transfers

The **P4-Adventure**'s seat moves forward 16 cm to ease transfers.

For slopes and descents

The **P4-Adventure** features a seat tilt to keep a horizontal position on steep slopes.



In your garden

Ideal for grass, gravel, stones, and all inaccessible paths to other wheelchairs.



In the city



The P4 definitively solves daily obstacle-crossing issues.

City travel is often difficult due to poor accessibility, particularly curbs. Despite recent improvements, many obstacles still block traditional wheelchairs. Running an errand can quickly become exhausting!

Their ability to climb obstacles higher than 12 cm lets you face anything with confidence!

Obstacle-climbing ability combined with an excellent turning radius make the **P4** surprisingly agile and practical in cities, where they gain a new vocation.

Climbing or descending a step is no longer an issue, and moving across uneven ground like cobblestones becomes easy.

Foldable roll bar (*)

Tilts forward or backward for easier vehicle access.





Road lighting + orange beacon (*)

Mounted on top of the roll bar, makes the wheelchair visible from far away in the dark.

Headrest

Integrated comfort headrest.

Seat with lateral support

Ensures trunk stability on rough paths.



Storage box (*)

Folds away with one arm to keep everything within reach.



Roll bar (*) and foldable backrest

For easier access into a vehicle.



Seat tilt (*)

Shifts center of gravity to secure steep descents.



Suspensions

Absorb vibrations and uneven ground.
Improve driving comfort and protect the back.



Cars fasteners

For car transport.

Luggage rack.

Low-pressure tires

Ensure optimal grip.



Anti-gutter armrests

Enhance arm stability on rocky
surfaces. Adjustable on 3 axes,
foldable backward.



R-Net electronics: 5 driving profiles.

Allows attendant control functions,
environmental control... Compatible with iOS –
Android + Bluetooth connection.

Forward seat adjustment (*)

Eases transfers, adjusts center
of gravity/gravité.



Controller mounted on swing away

With integrated cup
holder.



Transfer bars

Make getting in and
out easier.



Elevating leg rests.



Footplate

Tiltable from 90° to 180°.
Independently adjustable right/left.
Foldable backward.



Main power switch

Facilitates winter storage and
maintenance.











Curb descent wheels











Technical Specifications

P4 - 4x4

**The all paths wheelchair
where 4 motors power 4 wheels**

	P4 Crossover	P4 Adventure		P4 Crossover	P4 Adventure
 Range (depending on use)	>25 km	>45 km	 Obstacle clearance	>12 cm	>12 cm
 Top speed	6 ou 10 km/h	6 ou 10 km/h	 Turning radius	1,21 m	0,75 m
 Climbing capacity	20%	20%	 Motor power	4X350 Watts	4X350 Watts
 Descending capacity	20%	20%	 2 waterproof batteries	85 Amp (AGM+)	120 ou 180(*) Amp (LifePO4)
 Side slope capacity	15%	15%	 Charger: 15-25 Ah recommended IEB.	15 Amp	15 ou 25 Amp

Technical Specifications - Dimensions

	P4 Crossover	P4 Adventure		P4 Crossover	P4 Adventure
 Weight of the wheelchair without battery with battery	129 kg/ 215 kg	163 kg/ 189 kg	 Length of wheelchair - without footplates - with footplates raised - with footplates lowered	117 cm 125 cm 145 cm	117 cm 125 cm 145 cm
 Max user weight	120kg	120kg	 Width	69 cm	69 cm
 Seat tilt (*) - without top case - with top case	N/A	47° 15°	 Height - seat base - front seat with cushion - folded backrest (transport)	53 cm 65 cm N/A	53 cm 65 cm 99 cm
 Footplate angle	90 à 180°	90 à 180°	 Roll bar unfolded	N/A	150 cm

4POWER4

info@powerstand.eu

www.4power4-powerstand.com

Tested according to ISO 7176
standards.



Conforms to Annex 1 of
Directive 93/42/EEC.



Distributed by

(*) Optional depending on model and country.