



**URSUS ACTIVE™**  
active wheelchair



**URSUS ACTIVE™** is a special wheelchair designed for daily rehabilitation of children with locomotor dysfunctions. Furthermore, this wheelchair may be very useful concerning orthopaedic disorders and lower limbs relief needed after surgery.

The wheelchair is equipped with anti-tippers wheels that prevent the patients to tip over backwards. Optionally, Ursus **Active™** can also be equipped with push handles for a caregiver.

Tilt-in-space function allows you to smoothly adjust the tilt angle of the seat which helps you to relief hips by moving the weight of body from pelvis area to the back.

Thanks to the **ClickClack®** system we can quickly detach the seat and attach it to other products that are in our offer (**Ursus™**, **Ursus Home™**, **Ursus Active™**, **Quark™**).



## IT IS RECOMMENDED IN THE FOLLOWING CASES:

- Cerebral palsy – various types and forms.
- Brain and spinal cord injuries with lower limbs paralysis and weakness, balance and motor coordination disorders with different etiology.
- Demyelinating disease (e.g. multiple sclerosis – SM).
- Meningocele.
- Muscles diseases with weakness and paralysis (dystrophy, myasthenia gravis, myositis).
- Genetic disorders.
- Degeneration and neurological diseases.
- Other condition connected with muscles weakness, paralysis and musculoskeletal disorders.
- After head and spinal cord injuries and spine trauma with paralysis

## STANDARD EQUIPMENT:

Light aluminium frame | Large, inflatable rear wheels in sizes from 22" to 26" | equipped with quick-release couplings and strings | 2 independent rear brakes | Front shock absorbers | Seat unit with system ClickClack™ system for easy stick in/out | Fully reclining backrest to make the stroller as a bed | Adjustable backrest | Anatomical lateral supports | 5-points safety belts with height adjustment | Height adjustable backrest | Angle adjustable backrest | Reflective safety marks | Pneumatic tilt-in-space system TILT-IQ® | Height and angle adjustable footrest | 5 points feet stabilizing belts with heel support and foot | abduction-adduction adjustments | Angle adjustable footplate | Breathable upholstery | Soft front bar with reflective safety marks

*Our passion – your independence!*



### SPECIAL FEATURES:



#### Stable structure

The device ensures adequate stability and safety of Patient



#### Reflective elements

The product has reflective elements to ensure your safety.



#### Lightweight construction

The product is made of light metal alloys which positively affects the low weight of the device.

# 24

#### 24 months warranty

We offer 24 months warranty on the product. We also provide after sales service.



#### Easy to use

Intuitive way of adjustment and easy assembly.



#### Movable seat

You can turn the seat easily so that the child may be facing towards caregiver or movement direction.



#### Bucket Seat

It is possible to change the tilt angle of the seat.



Ergonomic push handle with height adjustment



Rear anti-tippers



Small swivel front casters helps to maneuver even on very small space

### STANDARD COLOURS:

#### BREATHABLE UPHOLSTERY:



Black  
Standard



Black & White  
Optional



Red  
Optional



Grass  
Optional



Sea  
Optional



Latte  
Optional

\*You can change upholstery colours between sizes without additional charges. Production time up to 14 days.

#### CITY UPHOLSTERY:



City Grass  
Size 1



City Rose  
Size 2



City Coffee  
Size 3



City Carbon  
Optional


*Our passion - your independence!*

Akces-MED sp. z o.o.  
Jasionka 955B  
36-002 JASIONKA

akces-med@akces-med.com  
tel. +48 (17) 864 04 70  
www.akces-med.com

Real colors may be different from the presented images.

version: 06/2016

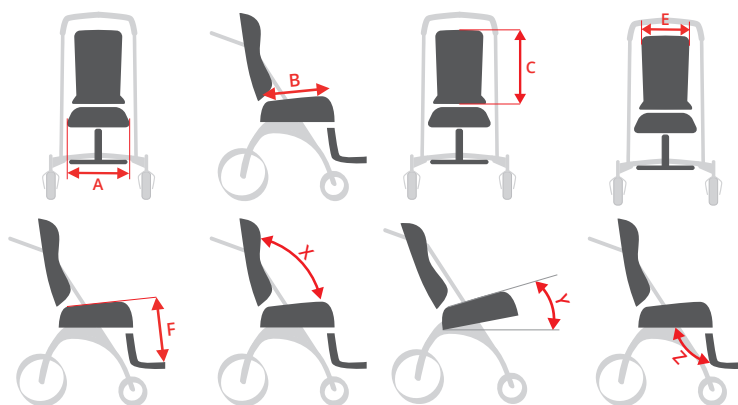

Symbol	Measurement	Unit	Size 1	Size 2	Size 3
A	Seat width	[cm]	19-25	22-30	24-34
B	Seat depth	[cm]	21-26	27-35	34-42
C	Backrest height	[cm]	44-55	49-59	55-67
E	Backrest width	[cm]	23	27	32
F	Footrests length	[cm]	13-29	16-40	16-43
	Footplate length	[cm]	20	20	20
	Push handle height	[cm]	85	85	85
X	Seat tilt angle (forward-facing)	[°]	90-180	90-180	90-180
Y	Seat tilt angle (rearward facing)	[°]	0-25	0-25	0-25
Z	Footrests tilt angle	[°]	90-180	90-180	90-180
	Maximum user weight	[kg]	35	35	50

Weight and dimensions of the device		Size 1	Size 2	Size 3
Width	[cm]	62	66	70
Length	[cm]	100	100	100
Height	[cm]	87	92	87
Weight	[kg]	18,5	20	22,5

Note: Measurement tolerance range is +/- 1 cm / ½ in.



### Safety



#### Soft cover for tray

Custom made item. production time 14 days.




#### Led Lights ( 2 pcs.)

Front and rear LED lights for increased safety.




### Positioning

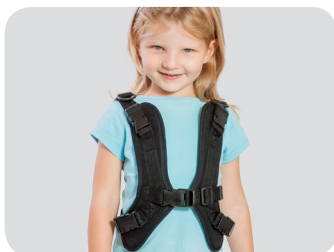

#### H vest

It acts as a trunk stabilizer. HOLD&PULL© system makes it easier to adjust the vest to the patient body. You can quickly detach the vest if necessary

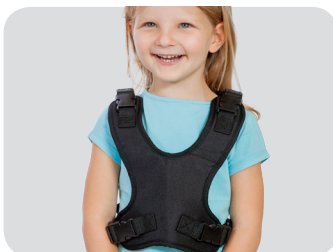

#### 4 points safety vest

It acts as a trunk stabilizer. Wide range of adjustment helps to fit it to the user's needs. Equipped with HOLD&PULL© system



#### 6 points safety vest

Wide range adjustable vest that acts as a trunk stabilizer. HOLD&PULL© system makes it easier to adjust the vest to the patient body.



#### Tight abduction belts

They can be easily attached to the chair by use of Velcro straps. Their use allows to keep legs in abduction and also prevents the patient from slipping off of the device.



#### Abduction-adduction thigh belts



#### Pelvic belt



#### Firm headrest made of foam



#### Head-neck support



#### Loins support



#### Trunk stabilizing belt

It acts as an additional trunk protection which helps to keep a proper sitting position. The belt is fastened by use of Velcro straps and its size can be adjusted to the user.



#### Trunk belt

Universal solution applicable to many of our products. It secures the user from falling out of the device. Anatomical shape and Velcro straps allow for a comfortable fit.



#### Abduction block

*Our passion - your independence!*

Real colors may be different from the presented images.

Akces-MED sp. z o.o.  
Jasionka 955B  
36-002 JASIONKA

akces-med@akces-med.com  
tel. +48 (17) 864 04 70  
www.akces-med.com

version: 06/2016



**Chest support**




**Narrowing pads (2 pcs) 6 cm**



**Narrowing pads (2 pcs) 10 cm**



## Comfort

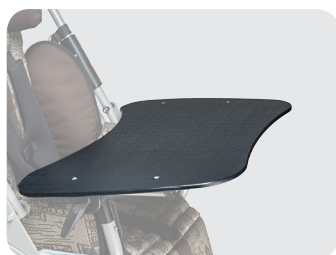

### Tray PU

It supports the therapy and allows you to play, eat, and work.



### Tray Plexiglas®

It supports the therapy and allows you to play, eat, and work.



### Headrest Cotton Cover

It helps to keep the headrest clean. Easy to assemble/disassemble, it can be washed.



### Cup Holder

An universal holder mounted to the frame. It enables easy and quick access to a bottle of drink.



**Decorated spokes guards (2 pcs)**



## Storage



**Lady's bag**



**Men's bag**



**Bag Delux**



### Oxygen tank holder

Oxygen tank case can be attached to the stroller by use of Velcro straps with adjustable length. Made to measure to perfectly fit the patient's oxygen tank.



### Build-up



Pneumatic Frame Quark™



Hi-Low Frame



Outdoor Frame



**Under seat storage basket**  
enables transporting small items in the chair. It is attached to the base of the chair.