



Care in motion

Automated pressure redistribution

Enhancing the users well-being

Supporting improved clinical outcomes

Toto[®] Cradle is an automated lateral turning system enabling prescribed lateral turning intervals, even when users are sleeping

With a powered control unit and a turning platform fitted beneath the users existing mattress, Toto Cradle automatically turns users at clinician-defined intervals from 30 minutes up to 4 hours, relieving pressure over their most vulnerable areas.

The new Toto Cradle is designed to offer enhanced comfort, support, and increased concordance

This system features a secondary tilt on the opposite side of the primary tilt, which work in parallel to cradle the user, providing enhanced comfort and support while minimising the risk of user migration across the support surface.

Launching
Toto Cradle
featuring
the

new



Cradle
Support
System

toto[®]
cradle



New Toto Cradle



Automated pressure redistribution

Toto Cradle provides effective pressure redistribution aiding pressure ulcer prevention and management.



Enhancing the users well-being

Toto Cradle is design to aid user well-being through supporting pressure area care management, preserving dignity, reducing pain experience compared to manual reposition, supporting end-of-life and palliative care, as well as greatly improving overall comfort and sleep quality.



Supporting improved clinical outcomes

Toto Cradle provides healthcare staff with peace of mind that the user is repositioned effectively at prescribed intervals whilst increasing their safety.


Designed to offer enhanced comfort, support, and increased concordance



Secondary Tilt*



Primary Tilt*

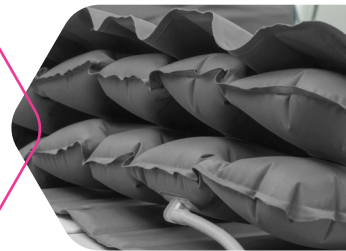

Clinically Effective Repositioning

Consistent and Controlled Turning

Promoting Safety

Improving Concordance

Adaptive Profiling Platform for all bed types

Easy to Use and Quick to Install

Inflation time of less than 3 minutes



Toto Touch Control Unit is compatible with Toto and Toto Cradle platforms

CPR Valve & Pipes

Toto Turning Platform & Cover

Toto Touch Control Unit

*Primary and secondary tilt functionality available on both side selections

Description	Dimensions	Weight	Order Code
Toto Cradle Platform	L 1850mm x W 745mm	6.9kg	5100000
Toto Touch Control Unit	L 240mm x W 100mm x H 200mm	2.2kg	4000000

For more information:

Call +44 (0)330 460 6030

Email: info@frontier-group.co.uk

Website: www.frontier-group.co.uk



Frontier Medical Group, Innova One, Tredegar Business Park, Tredegar, Blaenau Gwent NP22 3EL United Kingdom.

CE © Copyright Frontier Therapeutics Limited, 2023 All rights reserved. GP0133_B82_01



Toto Cradle lateral patient turning system has a two-year warranty



250kg User Weight Limit

toto
cradle