



**Yes, you can.®**

## **Alber viaplus**

### **Pushing comfort chairs the easy way**

Comfort wheelchairs offer users great comfort but sometimes pushing these wheelchair can be a strenuous job for the attendant. **Alber viaplus**, the push and brake aid for Invacare® comfort wheelchairs, offers attendants the assistance they need when pushing this type of wheelchair. What's more, the automatic braking function eliminates downhill worries and enables the safe negotiation of slopes and ramps.





### ➤ Fits perfectly to *Rea™ Azalea®* and *Clematis®*

Unlike many other pushing aids, the viaplus does not impede the adjustability of a tilt in space wheelchair (with a seat height of 45 cm). All of the wheelchair's functions remain accessible. That's why the viaplus can simply be left on the wheelchair after a walk without restricting the use of its functions.

### ➤ Easy to operate

The handy controls of the **Alber viaplus** fits onto the handle of the wheelchair without impeding its functions. The wheelchair is pushed forwards by light thumb pressure on the operating lever, making fatigue-free operation possible even on longer journeys. The chair is put into reverse by pulling the operating lever upwards.

### ➤ One-hand control

Controlling the **viaplus** is as easy as operating a single, handy rotary switch. This switches the drive on or off, sets the required speed or activates the freewheel mode. viaplus is therefore exceptionally easy to operate for older users.

### ➤ Personal freedom

**viaplus** not only benefits the attendant: It also gives the wheelchair user the opportunity for a fuller social life because frequent excursions with the wheelchair are no longer a problem now that handling the wheelchair is so much easier for the attendant. **viaplus** is also used in residential homes. With it, pushing wheelchairs becomes an effortless task and lightens the load for the staff.



## Features



### Drive unit

With puncture proof solid rubber tyres, automatic spring system for constant traction, pedal for lowering/lifting the wheels.



### Control unit

Central, handy rotary switch for simple control of all functions.



### Battery pack

Powerful battery pack with charging socket.



### Battery charger

Automatic charger with indicator light for charging status.



### Bracket set for wheelchair

Available for Rea™ **Clematis®** and **Azalea®**.



### Versatility

No limitation of seat/backrest adjustability (seat height = 45 cm).


## Options



### Battery box

Solid battery box for carrying the battery pack.

## Technical data

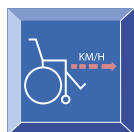

max. 5,5 km/h



16 km




24V/12Ah



9 kg



10%

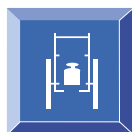

20 kg  
(drive unit: 10,7 kg)



**Azalea®**: 135 kg,  
**Clematis®**: 125 kg

All information quoted is believed to be correct at time of print. Invacare® reserves the right to alter product specifications without prior consultation. Ulrich Alber GmbH is a company and a brand of Invacare®.

Viaplus- EU - 06/10

© Invacare® International Sàrl 2010. All Rights Reserved.



Ulrich Alber GmbH  
Vor dem Weißen Stein 21 · 72461 Albstadt · Germany  
Tel.: +49 7432 2006 0 · Fax: +49 7432 2006 298  
Email: [info@alber.de](mailto:info@alber.de) · [www.alber.de](http://www.alber.de)

