

■ Head protection series Starlight®





Starlight®

Here you can see the range of atoform products, Europe's largest manufacturer of head protection for persons with special needs. Each product clearly displays the dedication, necessary experience and the high craftsmanship displayed during production, to achieve the most comfortable fit.

Meticulous choice of materials is taken for granted and well thought out production methods our demand.

The dialog with you is the foundation of sensible innovations.

Are you looking for something out of the ordinary?

Improvements and individual solutions are targeted and reached. We are grateful for your suggestions and wishes.

Our combined effort helps to help!

Your ATO FORM Team



Article number and
prices see
pricelist/order form

Starlight® Secure Evo



Elegance – with security!

The Starlight Secure Evo is based on the experiences and knowledge of the Starlight Secure Leather.

This model is distinguished by its improved materials for an even higher wearing comfort. High quality workmanship and elaborate workmanship, with much loving attention to detail are visible here. The upper leather is high-quality, toxin-free coloured full grain leather.

The closed upper side offers additional protection for those at risk of falling. Optional the head protection can be individually colour designed or can be offered with mixed materials.

With use of the Velcro sunshield the head protection becomes the "sporty" variant.

Accessories:

- Chin protection
- Forehead protection
- Neck protection
- Velcro fastener
- Fix lock fastener
- Strengthened thorn buckle



*Starlight® Secure Evo
with forehead and
neck protection
in various colour
(option)*



*Starlight® Secure Evo
with Velcro sunshield
and integrated denim
(option)*



Technical details: Starlight® Secure Evo

Size	1	2	3	4	5
Head circumference	47 – 49 cm	50 – 52 cm	53 – 55 cm	56 – 59 cm	60 – 62 cm
Padding	soft/firm foam combination				
Weight	approx. 240 g	approx. 260 g	approx. 270 g	approx. 290 g	approx. 310 g

Head protection



Starlight® Secure Leather



Elegance – with security!

The entire design is identical to the Starlight Secure head protection and is fitted with high-quality full grain cowhide.

The leather is vegetable tanned, free of chrome tanning agents and toxic substances, according to the PCP – prohibition laws.

The colouring is obtained by natural materials and processes, without use of AZO colours.

The inner lining consists of roughened calf leather. It is breathable, heat and moisture regulating, very kind to the skin, mineral nappa soft tanning in a barrel using high quality dyes.

Accessories:

- Chin protection
- Forehead protection
- Neck protection
- Fix lock fastener
- Strengthened thorn buckle

*Starlight® Secure Leather
with chin protection, forehead protection
and neck protection (Accessories)*



**Further colour requirements?
Certainly!**



Technical details: Starlight® Secure Leather

Size	1	2	3	4	5
Head circumference	47 – 49 cm	50 – 52 cm	53 – 55 cm	56 – 59 cm	60 – 62 cm
Padding	soft/firm foam combination				
Weight	approx. 150 g	approx. 180 g	approx. 180 g	approx. 200 g	approx. 240 g

Starlight® Easy

Especially light weight protector for a high definition fit. It adjusts any head shape, due to the special high-tech foam. And thus protects the cranium through well thought out anatomical adjustability.

The workmanship meets highest quality requirements. The protector is offered two toned, with free choice of colour.

Simple and brilliant adjustment system for the ear straps. Thought out excellence of quality.

Accessories:

- Chin protection
- Contrasting seams
- Padded chin strap



Elegant proportions and safe



Ear line adjustable



Technical details:

	Starlight® Easy
Forehead protection thickness	1,5 cm
Head circumference	according to size
Padding	hard foam
Outer material	leather
Lining	leather
Weight depends on size	approx. 120 – 260 g



When ordering please use the ATO FORM sheet for needs assessment "Body protection"!

You will find it in our price list or on our website www.ato-form.com.

Head protection



Starlight® Varia



The best head protection is the one which is gladly worn. Combined colours, special fabrics and sizes are available on request, for a small additional surcharge.

Growth-related or after surgery, it is possible that the shape of the head and therefore certain measurements (such as the head circumference) change in a short period of time.

Through its ingeniously construction the head protection Starlight Varia is individually adjustable in head circumference and straps making it ideal to adapt to the individual measures of the person in need of care.

The protective effect in the occipital area is achieved by the individual padding chambers made out of shock absorbing high-tech foam.

Through the adjustable possibilities, the head protection is strengthened in the occipital area. By this, any impact forces when falling backwards are spread to the entire construction of the occipital area.

The Starlight® Varia is made out of high quality toxin-free coloured full-grain leather.

Fasteners:

- Velcro fastener
- Fix lock fastener



**Further colour requirements?
Certainly!**



Technical details: Starlight® Varia

Size	S	M	L	XL
Forehead protection thickness	approx. 10 mm			
Head circumference	48 – 51 cm	52 – 55 cm	56 – 59 cm	60 – 63 cm
Padding	High-tech foam			
Outer material	leather			
Lining	leather			
Weight	approx. 310 g	approx. 330 g	approx. 350 g	approx. 370 g

Starlight® Flex



You expect precision? Okay! We are flexible...

Due to illness or growth-related, after surgery – for various reasons – the shape of the head can change.

Owing to its ingenious construction, the head protection Starlight® Flex is individually adjustable in many areas, making it ideal to adapt to the individual measures of the user.

Your customer deserves a risk-free Starlight® Flex without metal parts or similar hazards.

Accessories:

- Chin protection
- Forehead protection
- Velcro sunshield
- Closed top
- Ear protection



Starlight® Flex
basic model



Black

Blue

Brown

Green

Red



Closed top, forehead,
chin- and ear protection



Alternative upper side and
velcro sunshield

Technical details: Starlight® Flex

Size	S	M	L	XL
Forehead protection thickness (without Accessories)	approx. 23 mm			
Head circumference	48 – 52 cm	52 – 56 cm	56 – 60 cm	60 – 63 cm
Padding	High-tech foam			
Outer material	leather			
Lining	spacer fabrics			
Weight	approx. 340 g	approx. 360 g	approx. 380 g	approx. 400 g

**Technical details:**

Forehead protection thickness	Starlight® Protect 1,5 cm
Circumference	according to size
Padding	hard foam
Outer material	leather
Lining	leather
Weight	approx. 260 g



Colour variants

We offer many extras as a standard, for example the following colours:



Black

Blue

Brown

Green

Red

Starlight® Protect

Greatly improved further advancement of the model „Orient“

The construction of the models Protect and Protect Plus is based on the sandwich principle, in order to absorb and spread any occurring hit and pushforces.

A hard foam shell is worked in to the protector on the outside of the occipital area. This is hit and splinter resistant, to absorb and spread any occurring forces in case of a backward impact.

The Protect, as well as Protect Plus, are made of high quality full-grain leather, and are fastened securely and firmly, by chin and neck straps. The fastening can be chosen, depending on the indication.

The protective effect in the cranium area is accomplished by longitudinal braces (open version), or through the closed perforated cap, made of shockabsorbing high-tech foam.

Does your idea of a head protection look different? You are welcome to send us your idea on the measurement sheet, drawn to scale. Or in a personal interchange.

Accessories for Starlight® Protect and Protect Plus:

- Chin, ear, cheek protection
- Closed top
- Fix lock, Velcro® fastener
- Magnet lock including key
- Strengthened thorn buckle
- Visor or Xtra face protection (only for Protect Plus)



Starlight® Protect Plus



An asset for safety and Quality!

The versions Protect and Protect Plus differ mainly in the size of the forehead protection: Protect 1,5 cm, Protect Plus 4 cm.

Choose between different shades of leather – Our high standards of quality apply also to coloured leather.

The calf leather used for the lining is perforated to make it permeable. This enables the exchange of moisture and heat, and thus makes it very friendly to the skin.



You will find our ATO FORM „Assessment of demand“ sheet in our price list and on our homepage www.ato-form.com as well.



*Xtra face protection
(only for Protect Plus)*



Chin protection



Ear and cheek protection



Closed top



Visor (only for Protect Plus)

Technical details:	Starlight® Protect Plus
Forehead protection thickness	4 cm
Circumference	according to size
Padding	hard foam
Outer material	leather
Lining	leather
Weight	approx. 280 g



Starlight® Protect Plus

with „Xtra“ face protection

Is even more security possible?

The Starlight® Protect Plus with Xtra face protection was specially developed for patients with a great need for face protection.

The „Xtra“ offers even more safety at increased risk of injury from frontal crashes or auto aggressive behaviour.

Particularly important is the use of non-perforated calf leather, as it provides additional breathability, warmth and moisture regulation in combination with skin friendliness.



Can't find the solution you're looking for here?

We invite you to visit our website where you will find more information about our extensive range of special designs.

We will also gladly assist you by phone:

Please phone 0049 6093 944-0

Starlight® Go Sport & Go Sport Plus

Best protection, elegantly implemented!

The construction of both models Go Sport and Go Sport Plus are based on the sandwich principle, for absorption and distribution of all shocks and impacts.

A hard foam shell is worked into the protector on the outside of the occipital area. This is hit and splinter resistant, to absorb and spread any occurring forces in case of a backward impact.

The Go Sport, as well as Go Sport Plus are made of high quality full-grain leather and sit firmly and secure on the head. Various fasteners and much more selectable, depending upon indication.

Choose between various shades of leather – Our high standards of quality apply also to coloured leather!

The calf leather used for the lining is perforated, which is breathable, moisture and heat regulating as well as being very skin friendly.

You will find our ATO FORM „Assessment of demand“ sheet in our price list and on our homepage www.ato-form.com as well.

Accessories:

- Closed top
- Fix lock
- Velcro fastener
- Magnet lock incl. key
- Strengthened thorn buckle




Technical details: Starlight®	Go Sport / Go Sport Plus
Forehead protection thickness	1,5 cm / 3,5 cm
Head circumference	according to size
Padding	hard foam
Outer material	leather
Lining	leather
Weight	approx. 250 g / 300 g

Quality in detail – the varie



ty as required • individuality to taste

ety of **ato** ...lust for life® **form** gmbh

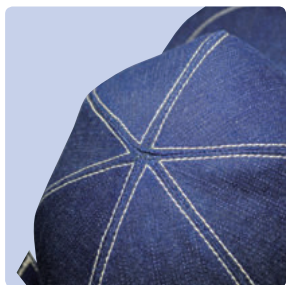

We want you to be in the best of hands. So we use: ■ Very soft yarn for highest wearing comfort ■ Solvent free glue ■ toxproof tanned full-grain leather ■ slightly perforated and roughened leather lining for an excellent microclimate ■ special highest quality foam ■ nickel free fasteners

Development, manufacturing and quality control carried out at our company headquarters in Germany.

Our in-house production gives us the possibility of **express production** - when the best protection is required urgently.



Starlight® Standard



This light weight head protection (about 100 g) has good deceleration properties.

The cotton jersey lining ensures good moisture circulation. It has a robust cover.

The adjustable chin strap offers excellent hold.

Starlight Standard is available in blue, red and brown. Each protector is manufactured individually and inch perfect.

For head circumference 45 to 62 cm available at standard price.

Washable at 40 Degrees.

Accessories:

- Chin protection
- Fix lock fastener
- Velcro fastener



Technical details:

	Starlight® Standard
Padding thickness	1,5 cm
Circumference	according to size
Padding	hard foam
Weight	approx. 100 g



Starlight Standard „Sport“ and Starlight Standard „Sport“ with closed top



Starlight Standard colour red, brown and blue

Starlight® Base



ATO FORM offers the most secure basis.

This head protection model is carefully made of high quality hardwearing cotton.

In addition to the lateral protection and the excellent fit, it offers excellent air circulation. Optionally, it can be supplemented with forehead and/or neck cushions, depending on the indications.

The relative protection in the cranium area is given through a crosscushion, made from shock absorbing foam material, that can be attached variably through Velcro fastener.


Adjustable chin strap and fasteners, which can be chosen according to indication, enable an excellent individual fit.

Accessories:

- Chin Protection
- Neck Protection
- Forehead Protection
- Fix lock fastener



Starlight® Base with chin guard, forehead protection and neck protection (Accessories)

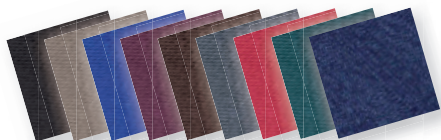

The best head protection is the one gladly worn. Colour combinations / special material / special sizes are available for all models on request, subject to a minimal surcharge.



Colour variants for Starlight® Base

Technical details: Starlight® Base

Size	1	2	3	4	5
Circumference	47 – 49 cm	50 – 52 cm	53 – 55 cm	56 – 59 cm	60 – 62 cm
Padding	soft/firm foam combination				
Weight	approx. 100 g	approx. 110 g	approx. 120 g	approx. 130 g	approx. 140 g



Starlight® Generation

Impressively lightweight protection

The Starlight® Generation models are padded with high-efficiency polyurethane foam. The reduced volume that characterises this head protection is made possible by the materials used. The protective effect is fully preserved.


Thus the Starlight® Generation models look more like “normal headwear” (e.g. cap, hat) and avoid the often stigmatising appearance of head protection.

The selected materials and cut guarantee optimal fit for lightweight protection.

The Starlight® Generation models can be optionally fitted with a fixed or removable chinstrap.

Accessories:

- Fixed or removable chinstrap, transparent dependant on model

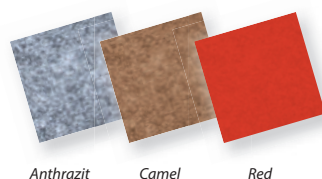

Technical details:

Model	Felt cap
Padding thickness	approx. 10 mm
Head circumference	according to size
Padding	Hightech-foam
Outer material	Felt
Lining	Cotton blend
Weight	ca. 160 g



It is possible to jazz up the Starlight® Generation Felt cap model with various accessories.


We are happy to advise you by telephone on 0049 6093 944-0 or through one of our rehab consultants.



Starlight® Generation – Colour contrasted seam stitching in denim look.



Starlight® Generation – Width adjustment

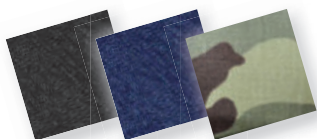

Jeans black Jeans blue Camouflage

Technical details:

Model	Baseball cap
Padding thickness	approx. 10 mm
Head circumference	according to size
Padding	Hightech-foam
Outer material	Cotton, mixed fabric
Lining	Cotton blend
Weight	ca. 160 g



Starlight® Baseball Cap – in special colour red with chinstrap or in denim look

Head protection



Starlight® Secure

*Starlight® Secure
with chin protection,
forehead protection
and neck protection
(accessories)*



The anatomical form of the “protector” with integrated protection of the temple and ear area, using high quality foam, ensures the excellent protection of the head. Hit and shock forces are absorbed through define pressure deformation.

The protection of the upper head area is accomplished through individually adjustable trapeze cushions, made of shock absorbing foam.


The head protection can be supplemented by forehead and/or neck cushions, which are fitted with Velcro strips. This ensures extraordinary protection. Individual and fitting is accomplished through adjustable chin straps as well as fasteners that can be individually chosen. Secure is covered by high quality cotton.

Accessories:

- Chin protection
- Forehead protection
- Neck protection
- Fix lock fastener



Colour variants

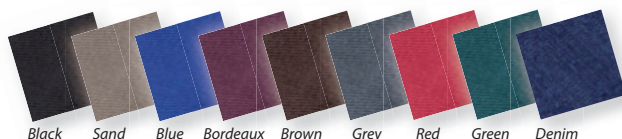

Require further accessories?

Certainly!

Technical details: Starlight® Secure

Size	1	2	3	4	5
Head circumference	47 – 49 cm	50 – 52 cm	53 – 55 cm	56 – 59 cm	60 – 62 cm
Padding	soft/firm foam combination				
Weight	approx. 100 g	approx. 100 g	approx. 110 g	approx. 120 g	approx. 130 g

Interesting aspects about the head protection program

You need a head protection for your child or another member of your family?

- The well thought out head protection line Starlight from ATO FORM can be acquired in your health product shop or orthopaedic shop. All models offer defined protection and high comfort levels.
- Model and material can be chosen according to indication. The different models are made of different special foams, chosen for specific requirements. You can also choose between different fasteners and colours for the head protection.
- In your health product - or orthopaedic shop, detailed specifications will be registered, using the ATO FORM sheet for needs assessment and indication. This is to ensure the excellent fit of the protector, to ensure its fulfilling its purpose to the full. This ensures, that the needs of the patient are not neglected, but are in first place.



Dimension sheet for Head protection Starlight® atoform

Comments (optional):

Frequency of falls: ☐ Often ☐ Often ☐ Very frequent

Fall direction: ☐ On front ☐ On back ☐ On side

Further notes: ☐ In case of falling injury to the head/neck

Therapist's name: _____ **Level of activity:** _____

General measurements:

Fig. 1

A: _____ cm (head circumference at forehead level, for the largest circumference)

B: _____ cm (distance from forehead to the top of the head)

C: _____ cm (distance from forehead to the top of the head, upper and lower part of the head)

D: _____ cm (head length, measuring and marking at the top of the head)

Position of the ears on the head:

Fig. 2

☐ Normal ☐ High ☐ Low

Skull perspectives:

Fig. 3

Fig. 4

Fig. 5

Fig. 6

Fig. 7

Fig. 8

Fig. 9

Fig. 10

Fig. 11

Fig. 12

Fig. 13

Fig. 14

Fig. 15

Fig. 16

Fig. 17

Fig. 18

Fig. 19

Fig. 20

Fig. 21

Fig. 22

Fig. 23

Fig. 24

Fig. 25

Fig. 26

Fig. 27

Fig. 28

Fig. 29

Fig. 30

Fig. 31

Fig. 32

Fig. 33

Fig. 34

Fig. 35

Fig. 36

Fig. 37

Fig. 38

Fig. 39

Fig. 40

Fig. 41

Fig. 42

Fig. 43

Fig. 44

Fig. 45

Fig. 46

Fig. 47

Fig. 48

Fig. 49

Fig. 50

Fig. 51

Fig. 52

Fig. 53

Fig. 54

Fig. 55

Fig. 56

Fig. 57

Fig. 58

Fig. 59

Fig. 60

Fig. 61

Fig. 62

Fig. 63

Fig. 64

Fig. 65

Fig. 66

Fig. 67

Fig. 68

Fig. 69

Fig. 70

Fig. 71

Fig. 72

Fig. 73

Fig. 74

Fig. 75

Fig. 76

Fig. 77

Fig. 78

Fig. 79

Fig. 80

Fig. 81

Fig. 82

Fig. 83

Fig. 84

Fig. 85

Fig. 86

Fig. 87

Fig. 88

Fig. 89

Fig. 90

Fig. 91

Fig. 92

Fig. 93

Fig. 94

Fig. 95

Fig. 96

Fig. 97

Fig. 98

Fig. 99

Fig. 100

Fig. 101

Fig. 102

Fig. 103

Fig. 104

Fig. 105

Fig. 106

Fig. 107

Fig. 108

Fig. 109

Fig. 110

Fig. 111

Fig. 112

Fig. 113

Fig. 114

Fig. 115

Fig. 116

Fig. 117

Fig. 118

Fig. 119

Fig. 120

Fig. 121

Fig. 122

Fig. 123

Fig. 124

Fig. 125

Fig. 126

Fig. 127

Fig. 128

Fig. 129

Fig. 130

Fig. 131

Fig. 132

Fig. 133

Fig. 134

Fig. 135

Fig. 136

Fig. 137

Fig. 138

Fig. 139

Fig. 140

Fig. 141

Fig. 142

Fig. 143

Fig. 144

Fig. 145

Fig. 146

Fig. 147

Fig. 148

Fig. 149

Fig. 150

Fig. 151

Fig. 152

Fig. 153

Fig. 154

Fig. 155

Fig. 156

Fig. 157

Fig. 158

Fig. 159

Fig. 160

Fig. 161

Fig. 162

Fig. 163

Fig. 164

Fig. 165

Fig. 166

Fig. 167

Fig. 168

Fig. 169

Fig. 170

Fig. 171

Fig. 172

Fig. 173

Fig. 174

Fig. 175

Fig. 176

Fig. 177

Fig. 178

Fig. 179

Fig. 180

Fig. 181

Fig. 182

Fig. 183

Fig. 184

Fig. 185

Fig. 186

Fig. 187

Fig. 188

Fig. 189

Fig. 190

Fig. 191

Fig. 192

Fig. 193

Fig. 194

Fig. 195

Fig. 196

Fig. 197

Fig. 198

Fig. 199

Fig. 200

Fig. 201

Fig. 202

Fig. 203

Fig. 204

Fig. 205

Fig. 206

Fig. 207

Fig. 208

Fig. 209

Fig. 210

Fig. 211

Fig. 212

Fig. 213

Fig. 214

Fig. 215

Fig. 216

Fig. 217

Fig. 218

Fig. 219

Fig. 220

Fig. 221

Fig. 222

Fig. 223

Fig. 224

Fig. 225

Fig. 226

Fig. 227

Fig. 228

Fig. 229

Fig. 230

Fig. 231

Fig. 232

Fig. 233

Fig. 234

Fig. 235

Fig. 236

Fig. 237

Fig. 238

Fig. 239

Fig. 240

Fig. 241

Fig. 242

Fig. 243

Fig. 244

Fig. 245

Fig. 246

Fig. 247

Fig. 248

Fig. 249

Fig. 250

Fig. 251

Fig. 252

Fig. 253

Fig. 254

Fig. 255

Fig. 256

Fig. 257

Fig. 258

Fig. 259

Fig. 260

Fig. 261

Fig. 262

Fig. 263

Fig. 264

Fig. 265

Fig. 266

Fig. 267

Fig. 268

Fig. 269

Fig. 270

Fig. 271

Fig. 272

Fig. 273

Fig. 274

Fig. 275

Fig. 276

Fig. 277

Fig. 278

Fig. 279

Fig. 280

Fig. 281

Fig. 282

Fig. 283

Fig. 284

Fig. 285

Fig. 286

Fig. 287

Fig. 288

Fig. 289

Fig. 290

Fig. 291

Fig. 292

Fig. 293

Fig. 294

Fig. 295

Fig. 296

Fig. 297

Fig. 298

Fig. 299

Fig. 300

Fig. 301

Fig. 302

Fig. 303

Fig. 304

Fig. 305

Fig. 306

Fig. 307

Fig. 308

Fig. 309

Fig. 310

Fig. 311

Fig. 312

Fig. 313

Fig. 314

Fig. 315

Fig. 316

Fig. 317

Fig. 318

Fig. 319

Fig. 320

Fig. 321

Fig. 322

Fig. 323

Fig. 324

Fig. 325

Fig. 326

Fig. 327

Fig. 328

Fig. 329

Fig. 330

Fig. 331

Fig. 332

Fig. 333

Fig. 334

Fig. 335

Fig. 336

Fig. 337

Fig. 338

Fig. 339

Fig. 340

Fig. 341

Fig. 342

Fig. 343

Fig. 344

Fig. 345

Fig. 346

Fig. 347

Fig. 348

Fig. 349

Fig. 350

Fig. 351

Fig. 352

Fig. 353

Fig. 354

Fig. 355

Fig. 356

Fig. 357

Fig. 358

Fig. 359

Fig. 360

Fig. 361

Fig. 362

Fig. 363

Fig. 364

Fig. 365

Fig. 366

Fig. 367

Fig. 368

Fig. 369

Fig. 370

Fig. 371

Fig. 372

Fig. 373

Fig. 374

Fig. 375

Fig. 376

Fig. 377

Fig. 378

Fig. 379

Fig. 380

Fig. 381

Fig. 382

Fig. 383

Fig. 384

Fig. 385

Fig. 386

Fig. 387

Fig. 388

Fig. 389

Fig. 390

Fig. 391

Fig. 392

Fig. 393

Fig. 394

Fig. 395

Fig. 396

Fig. 397

Fig. 398

Fig. 399

Fig. 400

Fig. 401

Fig. 402

Fig. 403

Fig. 404

Fig. 405

Fig. 406

Fig. 407

Fig. 408

Fig. 409

Fig. 410

Fig. 411

Fig. 412

Fig. 413

Fig. 414

Fig. 415

Fig. 416

Fig. 417

Fig. 418

Fig. 419

Fig. 420

Fig. 421

Fig. 422

Fig. 423

Fig. 424

Fig. 425

Fig. 426

Fig. 427

Fig. 428

Fig. 429

Fig. 430

Fig. 431

Fig. 432

Fig. 433

Fig. 434

Fig. 435

Fig. 436

Fig. 437

Fig. 438

Fig. 439

Fig. 440

Fig. 441

Fig. 442

Fig. 443

Fig. 444

Fig. 445

Fig. 446

Fig. 447

Fig. 448

Fig. 449

Fig. 450

Fig. 451

Fig. 452

Fig. 453

Fig. 454

Fig. 455

Fig. 456

Fig. 457

Fig. 458

Fig. 459

Fig. 460

Fig. 461

Fig. 462

Fig. 463

Fig. 464

Fig. 465

Fig. 466

Fig. 467

Fig. 468

Fig. 469

Fig. 470

Fig. 471

Fig. 472

Fig. 473

Fig. 474

Fig. 475

Fig. 476

Fig. 477

Fig. 478

Fig. 479

Fig. 480

Fig. 481

Fig. 482

Fig. 483

Fig. 484

Fig. 485

Fig. 486

Fig. 487

Fig. 488

Fig. 489

Fig. 490

Fig. 491

Fig. 492

Fig. 493

Fig. 494

Fig. 495

Fig. 496

Fig. 497

Fig. 498

Fig. 499

Fig. 500

Fig. 501

Fig. 502

Fig. 503

Fig. 504

Fig. 505

Fig. 506

Fig. 507

Fig. 508

Fig. 509

Fig. 510

Fig. 511

Fig. 512

Fig. 513

Fig. 514

Fig. 515

Fig. 516

Fig. 517

Fig. 518

Fig. 519

Fig. 520

Fig. 521

Fig. 522

Fig. 523

Fig. 524

Fig. 525

Fig. 526

Fig. 527

Fig. 528

Fig. 529

Fig. 530

Fig. 531

Fig. 532

Fig. 533

Fig. 534

Fig. 535

Fig. 536

Fig. 537

Fig. 538

Fig. 539

Fig. 540

Fig. 541

Fig. 542

Fig. 543

Fig. 544

Fig. 545

Fig. 546

Fig. 547

Fig. 548

Fig. 549

Fig. 550

Fig. 551

Fig. 552

Fig. 553

Fig. 554

Fig. 555

Fig. 556

Fig. 557

Fig. 558

Fig. 559

Fig. 560

Fig. 561

Fig. 562

Fig. 563

Fig. 564

Fig. 565

Fig. 566

Fig. 567

Fig. 568

Fig. 569

Fig. 570

Fig. 571

Fig. 572

Fig. 573

Fig. 574

Fig. 575

Fig. 576

Fig. 577

Fig. 578

Fig. 579

Fig. 580

Fig. 581

Fig. 582

Fig. 583

Fig. 584

Fig. 585

Fig. 586

Fig. 587

Fig. 588

Fig. 589

Fig. 590

Fig. 591

Fig. 592

Fig. 593

Fig. 594

Fig. 595

Fig. 596

Fig. 597

Fig. 598

Fig. 599

Fig. 600

Fig. 601

Fig. 602

Fig. 603

Fig. 604

Fig. 605

Fig. 606

Fig. 607

Fig. 608

Fig. 609

Fig. 610

Fig. 611

Fig. 612

Fig. 613

Fig. 614

Fig. 615

Fig. 616

Fig. 617

Fig. 618

Fig. 619

Fig. 620

Fig. 621

Fig. 622

Fig. 623

Fig. 624

Fig. 625

Fig. 626

Fig. 627

Fig. 628

Fig. 629

Fig. 630

Fig. 631

Fig. 632

Fig. 633

Fig. 634

Fig. 635

Fig. 636

Fig. 637

Fig. 638

Fig. 639

Fig. 640

Fig. 641

Fig. 642

Fig. 643

Fig. 644

Fig. 645

Fig. 646

Fig. 647

Fig. 648

Fig. 649

Fig. 650

Fig. 651

Fig. 652

Fig. 653

Fig. 654

Fig. 655

Fig. 656

Fig. 657

Fig. 658

Fig. 659

Fig. 660

Fig. 661

Fig. 662

Fig. 663

Fig. 664

Fig. 665

Fig. 666

Fig. 667

Fig. 668

Fig. 669

Fig. 670

Fig. 671

Fig. 672

Fig. 673

Fig. 674

Fig. 675

Fig. 676

Fig. 677

Fig. 678

Fig. 679

Fig. 680

Fig. 681

Fig. 682

Fig. 683

Fig. 684

Fig. 685

Fig. 686

Fig. 687

Fig. 688

Fig. 689

Fig. 690

Fig. 691

Fig. 692

Fig. 693

Fig. 694

Fig. 695

Fig. 696

Fig. 697

Fig. 698

Fig. 699

Fig. 700

Fig. 701

Fig. 702

Fig. 703

Fig. 704

Fig. 705

Fig. 706

Fig. 707

Fig. 708

Fig. 709

Fig. 710

Fig. 711

Fig. 712

Fig. 713

Fig. 714

Fig. 715

Fig. 716

Fig. 717

Fig. 718

Fig. 719

Fig. 720

Fig. 721

Fig. 722

Fig. 723

Fig. 724

Fig. 725

Fig. 726

Fig. 727

Fig. 728

Fig. 729

Fig. 730

Fig. 731

Fig. 732

Fig. 733</

Head protection



Starlight® Aqua



Until now it was not possible to wear the head protection when bathing/showering, or for therapeutic purposes or water based leisure activities. Practice shows however, that also in these areas protection of the head is very important. Therefore the head protection Starlight® Aqua was solely designed for this area.

The padding of the Aqua consists of a special integral foam, which is soft and elastic and therefore shock absorbent. Thanks to the careful selection of materials a high wearing comfort is achieved. The many positive properties of the neoprene used (high dynamic load capacity, weather and ozone resistance, high elasticity and impact resilience as well as resilience against fungi and bacteria, quick-drying, waterproof, easy to clean and disinfect etc.) makes this head protection a valuable aid.

Through the anatomical shape the entire head is protected. The adjustability of the chin strap and upper head ensures that the head protection is neither too tense nor too loose.

The cleaning of the head protection is also no problem by hand washing up to 30°.



*Individuell verstellbare
Trapetzipolster*



Chin protection (accessory)

Accessory:

- Chin Protection

Fasteners:

- Velcro® fastener
- Fix lock fastener

Technical details: Starlight® Aqua

Size	1	2	3	4	5
Head circumference	47 – 49 cm	50 – 52 cm	53 – 55 cm	56 – 59 cm	60 – 62 cm
Padding	extreme soft integral-foam				
Outer material	neoprene, breathable and non-toxic				
Weight	approx. 130 g	approx. 150 g	approx. 170 g	approx. 190 g	approx. 210 g

Starlight® Snooze

An exceptional light protector with good deceleration qualities and a highly defined fit. Extra soft, breathable leather and soft yarns are used for processing.

The comfortable soft padding made out of a combination of various foam types, will guaranteed prevent unpleasant pressure whilst sleeping.

Technical details:	Starlight® Snooze
Padding thickness	4 cm
Circumference	according to size
Padding	hard foam + marshmallow foam + cotton wool
Outer material	leather
Lining	leather
Weight	approx. 270 g



Accessories:

- Chin protection
- Strengthened thorn buckle

Starlight® Active

The use of modern communication methods, or more traditional tools like paint brush or pencil, is made possible through a pressure rod, which is adjustable in length and tilt.

This light weight head protection is available in two sizes, and easy to adjust to differing head circumferences. It is made of soft leather (Weight approx. 250 g).

Product contents:

- Head protector with pressure rod

Indication:

- Well preserved motor activity of the head



Technical details: Starlight® Active

Size	1	2
Head circumference	up to 53 cm	from 53 cm

Head protection

Trial fitting Your head protection supply for the best fit

Dimension sheet for Head protection Starlight® **ato form**

Commission (anonymous) ☐ ☐ ☐

Frequency of falls ☐ infrequent ☐ often ☐ very frequent

Fall direction ☐ forward ☐ backward ☐ left ☐ right

Is there risk of head injury or self-harming behaviour? ☐ Yes ☐ No

Disease profile Age: _____ Level of activity: _____

General measurements Fig. 1

A _____ cm (Head circumference at eyebrow height, above the occipital protuberance)

B _____ cm (Distance between the eyebrow to upper edge of the occipital protuberance)

C _____ cm (Over the highest point of the head, beginning and ending at the line of head circumference measurement A)

D _____ cm (Under the chin, beginning and ending at the line of head circumference measurement A)

*Occipital protuberance

Position of the ears on the head Fig. 2

☐ normal ☐ low ☐ high

Skull perspectives Fig. 3

Back view:

Side view:

Front view:

The head dimensions are indicated to compare feasible helmet types according to the distribution tables.

Model incl. Additional equipment, Accessories: _____

Colour: _____

Pattern: _____

Attention

In the case of anatomical abnormalities, describe the shape of the head as precisely as possible. Photos are very helpful here.

Identify the type of the ears are located on the line of measurement A.

If the positioning of the ears differs from this, make certain that measurement C is taken across the top of the head, beginning and ending at the line of measurement A.

In the case of asymmetrical ear positions, please give a color and detailed description in the box labelled "Comments".

Comments: _____

If necessary, please print and attach the forehead graphic and shape in the adjacent illustration as applicable, note any other dimensions.

Address/Stamp or Customer number: _____


Name/Company: _____

Address: _____

Phone: _____

Customer number: _____

Designation: _____

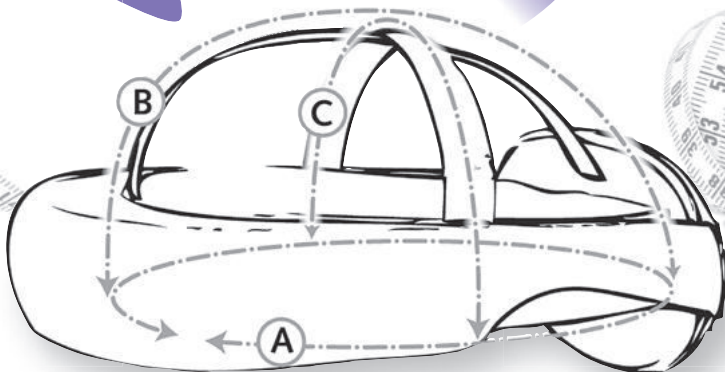

...enhanced wearing comfort and secure protection.

Make use of the opportunity, to test the fit of the head protection firstly as a trial fitting on your patient before completion of the head protection.

We will gladly send you the trial fitting to adapt at cost price, no matter which model.

Small corrections to the trial fitting are unproblematic possible. In difficult cases, you will also find the correct dimension.

An easy and secure control of the fit is offered.

We gladly will advise you if you have questions by telephone +49 6093 944-0.




To avoid errors, we advise to **only use our current dimension sheet**



Comfort and protection on highest level...



Theme brochures from ATO FORM

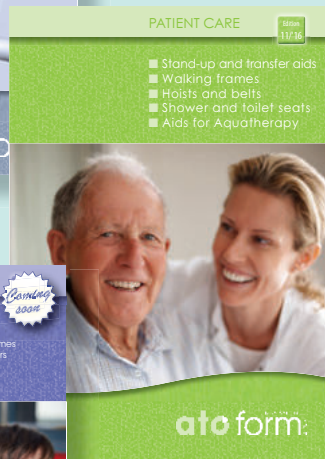

ato form ...just for life® gmbh

www.ato-form.com

Natasha Kilburn
Export

Phone: +496093 944 170
Mobil: +49151 14830125
export@ato-form.com

