

MagniLink Zip Premium

Versatile and easy to use!

The portable and foldable MagniLink Zip is well known for its image quality and user-friendly operation. Yet it provides all the functions offered in an advanced, stationary video magnifier.

Light and easy

MagniLink Zip is easy to fold up and down and can easily be stowed away or brought along on a trip. When combined with a battery it becomes a really neat travelling companion.

User-friendly and ergonomic

The control panel with tactile buttons and knobs is very easy to operate and the system also offers an advanced mode for advanced users.

The monitor is very flexible and can be set to the height and angle that suits you the best. By adjusting the monitor, the way that your head and eyes are in a 20-degree angle downwards, you will sit more ergonomically and reduce the tension on your neck and shoulders.


Your MagniLink Zip comes in a practical carrying bag that will ease transportation. The bag has orange details which makes it easier to find zippers and pockets.

Built-in distance camera

The camera can easily be folded up – and you have a distance camera in a second! The distance camera is, for example, perfect in a school environment when the student needs to see what the teacher is doing at the board. You can adjust the distance camera both horizontally and vertically. With the distance camera, you also have access to a really great mirror function that is perfect, when you for instance want to do makeup or trim your beard.




t always fascinated me. For
ars I worked hard to achieve
orts. I photographed everyt
ed and put lots of time in th
er. Sports took me out in th
orld and I was travelling to b

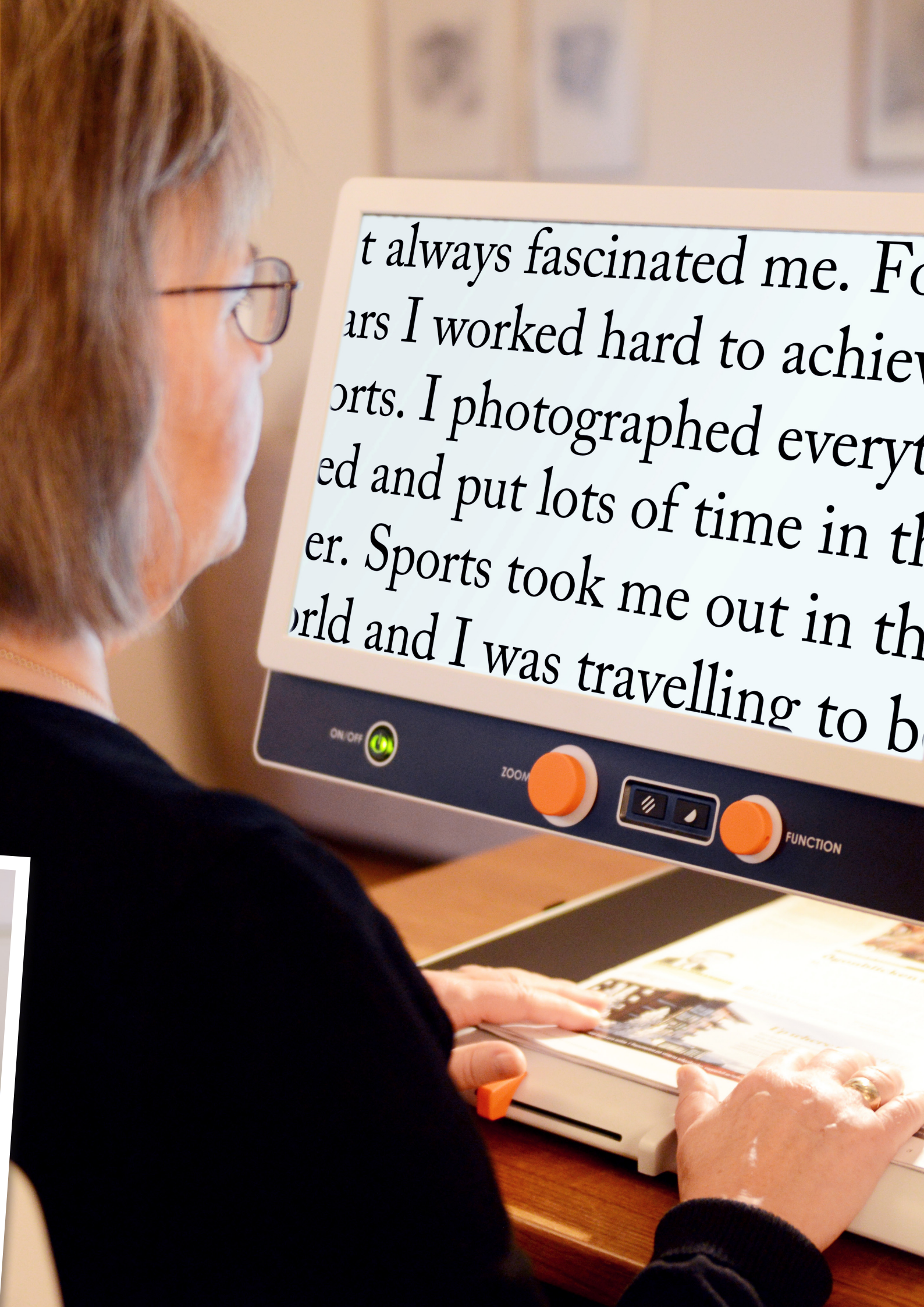

Make it yours!

Customize the MagniLink Zip for your specific needs. It is suitable for most tasks at work, in school and at home. The system is a great support when, for example, you want to read different documents, see the board in the classroom or do handicraft.

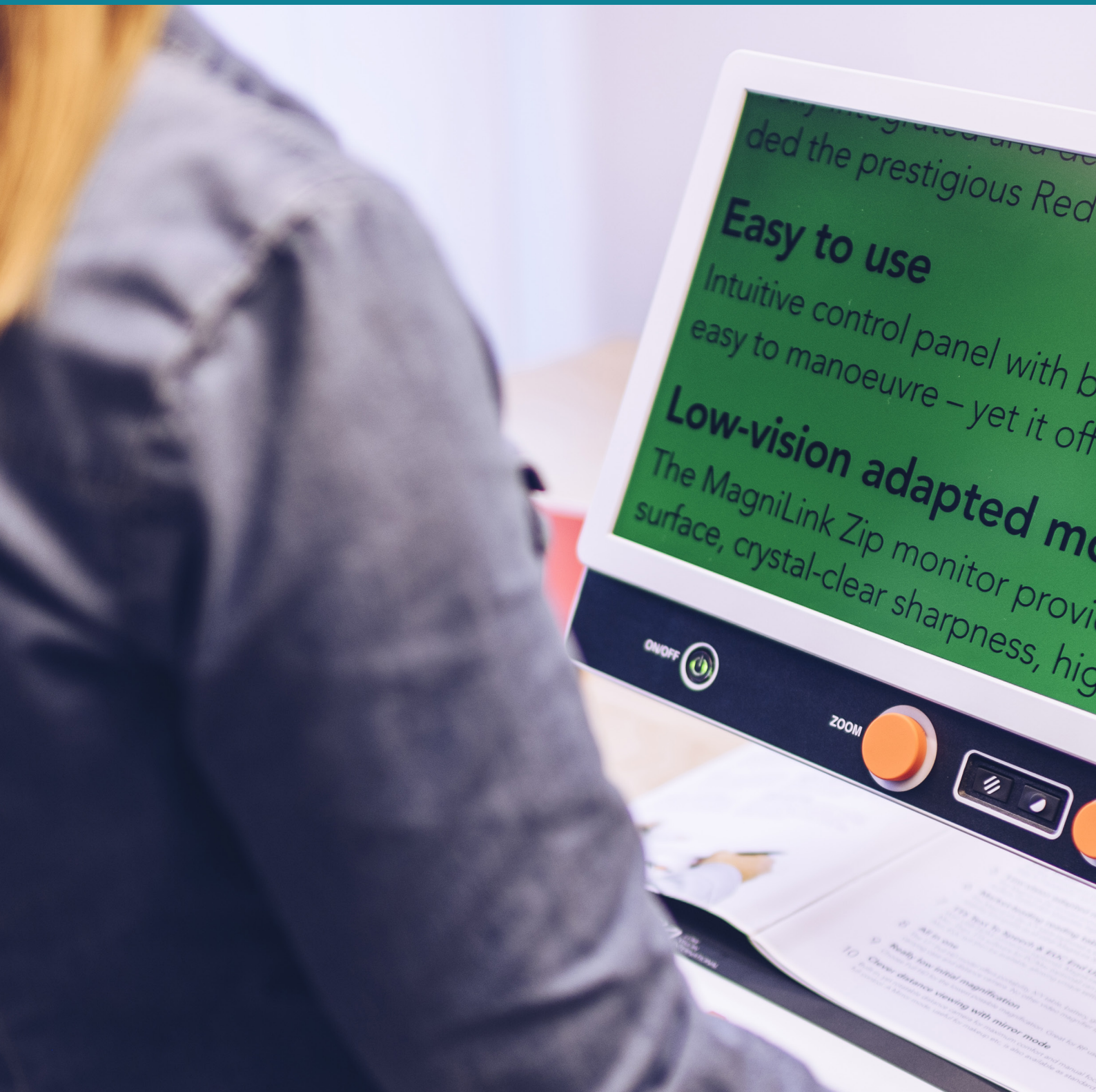



CHOOSE A HD OR
FULL HD CAMERA

CHOOSE
MONITOR SIZE

CHOOSE A UNIT WITH OR
WITHOUT READING TABLE



Outstanding image quality

MagniLink Zip is available with HD or Full HD camera. The HD camera gives you a good picture and is suitable for most purposes. Full HD provides the sharpest image quality ever and is, for instance, the ideal choice for RP users as it offers the best contrast and low initial magnification.

Low-vision adapted monitors

Choose between two screen sizes, 13" or 17", for your MagniLink Zip. The smaller screen is ideal if you want to be more flexible and need to replace and transport the unit quite frequently. Are you however looking for a more stationary system, but still wants to be able to fold and move the unit as smoothly as possible. Then the system with 17" screen and the integrated reading table is the obvious choice.

Convenient reading table

The integrated reading table is stable yet thin and has friction brakes. The table is a perfect tool when you move documents while reading, especially at high magnification. Observe that the integrated reading table is factory-installed and needs to be ordered from the start. If you order a unit without the integrated reading table, you can instead choose the smaller dockable table with friction brakes, MLRT-XY. This table can be ordered separately.



Computer connection

With USB/HDMI connection you can connect your MagniLink Zip to a computer (PC/Mac/Chromebook). Then you can view the camera image on your computer screen instead. One of the advantages of computer connection is that you can choose a text-to-speech function (TTS). To have text read aloud enables increased reading speed as well as reading endurance.

Choose a unit with battery

If you need to move around your MagniLink Zip often, or perhaps use it during travel, it may be smart to choose a unit with battery. Then you become more mobile and flexible in your use – perfect for people on the go!

Red Dot Award

MagniLink Zip is awarded the Red Dot Award for its user-friendly design. The reward shows that MagniLink Zip not only lives up to the high demands made by the users but also that the system is the best in its category.





Technical specifications

Model	MagniLink Zip 13"	MagniLink Zip 17"
Camera	HD: 720p, FHD: 1080p	HD: 720p, FHD: 1080p
Screen resolution	1600x900	1920x1080
Magnification	HD: 1.3-35x, FHD: 1.1-35x	HD: 1.7-45x, FHD: 1.4-45x
Battery, operation time	Up to 5 h	Up to 4.5 h
Battery, charging time	2.5 h	2.5 h
Dimensioner, folded	33x45x9 cm (13x17.7x3.5 in)	44x47x12 cm (17.3x18.5x4.7 in)
Dimensioner, unfolded	33x45x51 cm (13x17.7x20 in)	44x47x56 cm (17.3x18.5x22 in)
Weight	from 4.4 kg (9.7 lbs)	from 5.1 kg (11.2 lbs)
Power source	Battery or AC	Battery or AC

MagniLink Zip Premium

LVI Low Vision International is recognized as one of the world's leading manufacturers of low vision aids with the mission to make daily life easier for people with visual impairments.

We are certified to ISO 14001 & 9001



LVI Low Vision International
Verkstadsgratan 5
352 46 Växjö, SWEDEN
Tel: +46 (0)470 72 77 00
Fax: +46 (0)470 72 77 25
E-mail: info@lvi.se
WWW.LVI.SE